

## cFluor<sup>®</sup> B548 Anti-Human CD20 (2H7)

PRODUCT DETAILS	
<b>Catalog Number:</b>	R7-20201 (100 tests) R7-20202 (25 tests)
<b>Reactivity:</b>	Human, Baboon, Capuchin Monkey, Chimpanzee, Cynomolgus, Pigtailed Macaque, Rhesus, Squirrel Monkey
<b>Clone:</b>	2H7
<b>Format:</b>	cFluor <sup>®</sup> B548*
<b>Isotype:</b>	Mouse IgG2b, κ
<b>Test Dilution:</b>	5 μL / test
<b>Application:</b>	Flow cytometry
<b>Formulation:</b>	Phosphate-buffered saline, pH 7.2, containing 0.09% sodium azide and 0.2% BSA (BSA Country of Origin USA)
<b>Storage:</b>	2-8°C and protected from light. <b>Do not freeze</b>

### PRODUCT DESCRIPTION

The 2H7 monoclonal antibody binds to human CD20, a 33-37 kDa, four-transmembrane phosphoprotein. CD20 is present on all stages of B cell development except for early pro-B cells or plasma cells to play a role in B-cell activation, proliferation and differentiation<sup>1</sup>. In some studies, CD20 is thought to act as a Ca<sup>2+</sup> channel in the plasma membrane<sup>2</sup>. It is also present on a T cell subset at low levels as well as some follicular dendritic cells. The antibody was conjugated to a fluorophore and purified by affinity chromatography.


### RECOMMENDED USAGE

Each lot of this antibody is quality control tested using flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 μL per 1 million cells in a staining volume of 100 μL. If whole blood is analyzed, then use 5 μL per 100 μL. It is recommended that users titrate the antibody to obtain the optimal result for their specific application.

Please briefly centrifuge the reagent vial before use.


Use appropriate personal protective equipment per the product safety data sheet when using this product.

PRODUCT DATA

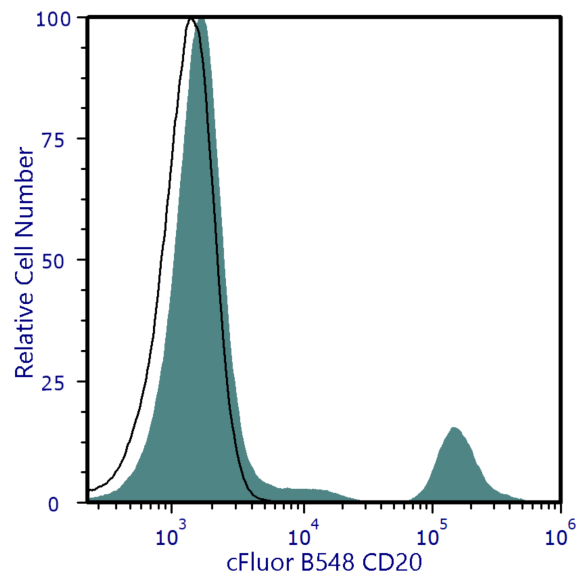

Human peripheral lymphocytes stained with cFluor® B548 CD20 (clone 2H7) (filled histogram) or mouse cFluor® B548 IgG2b, κ isotype control (open histogram).

Spectral signature of cFluor® B548 from a Cytek® Aurora 5 laser system equipped with 355, 405, 488, 561 and 640 nm lasers using CytekAssaySetting.

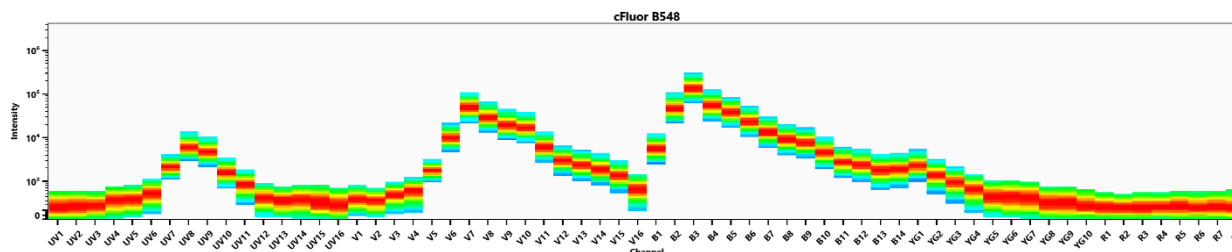

REFERENCES

1. Tedder TF, et al. 1994. Immunol Today. 15(9):450–4
2. Putney JW Jr. 1999. Proc Natl Acad Sci U S A. 96(26):14669-14671

**For Research Use Only. Not intended for use in diagnostic procedures.**

\*cFluor™ B548 is equivalent to CF® 514 manufactured and provided by Biotium, Inc. under an Agreement between Biotium and Cytek (LICENSEE). The manufacture, use, sale, offer for sale, or import of the product is covered by one or more of the patents or pending applications owned or licensed by Biotium. The purchase of this product includes a limited, non-transferable immunity from suit under the foregoing patent claims for using only this amount of product for the purchaser’s own internal research. No right under any other patent claim, no right to perform any patented method, and no right to perform commercial services of any kind, including without limitation reporting the results of purchaser’s activities for a fee or other commercial consideration, is conveyed expressly, by implication, or by estoppel.