

## cFluor<sup>®</sup> B548 Anti-Human CD45 (HI30)

PRODUCT DETAILS	
<b>Catalog Number:</b>	R7-20291 (100 tests) R7-20292 (25 tests)
<b>Reactivity:</b>	Human, Chimpanzee
<b>Clone:</b>	HI30
<b>Format:</b>	cFluor <sup>®</sup> B548*
<b>Isotype:</b>	Mouse IgG1, $\kappa$
<b>Test Dilution:</b>	5 $\mu$ L / test
<b>Application:</b>	Flow cytometry
<b>Formulation:</b>	Phosphate-buffered saline, pH 7.2, containing 0.09% sodium azide and 0.2% BSA (BSA Country of Origin USA)
<b>Storage:</b>	2-8°C and protected from light. <b>Do not freeze</b>

### PRODUCT DESCRIPTION

The HI30 monoclonal antibody binds to isoforms of human CD45, a 180-240-kDa type I membrane glycoprotein, also known as LCA (leukocyte common antigen). CD45 is the first and prototypic receptor-like protein tyrosine phosphatase that expresses on all human leukocytes<sup>1,2</sup>. It is absent on mature erythrocytes, platelets and non-hematopoietic cells. The antibody was conjugated to a fluorophore and purified by chromatography.


### RECOMMENDED USAGE

Each lot of this antibody is quality control tested using flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5  $\mu$ L per 1 million cells in a staining volume of 100  $\mu$ L. If whole blood is analyzed, then use 5  $\mu$ L per 100  $\mu$ L. It is recommended that users titrate the antibody to obtain the optimal result for their specific application.


Please briefly centrifuge the reagent vial before use.

Use appropriate personal protective equipment per the product safety data sheet when using this product.

PRODUCT DATA

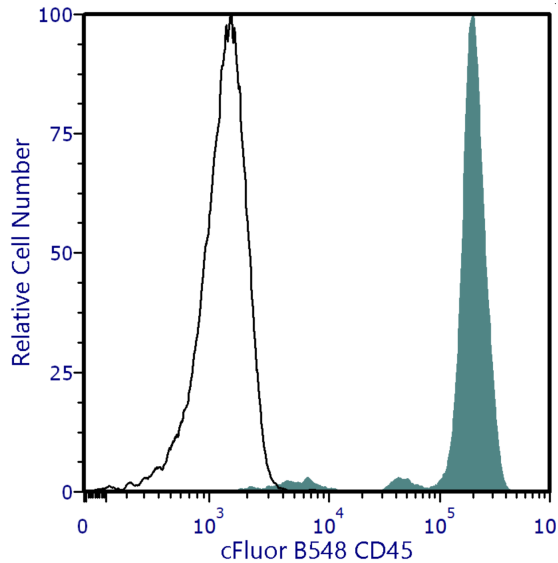

Human peripheral blood was stained with cFluor® B548 Anti-Human CD45 (clone 2D1) (filled histogram) or cFluor® B548 mouse IgG1, κ isotype control (open histogram). Data shown was gated on lymphocytes.

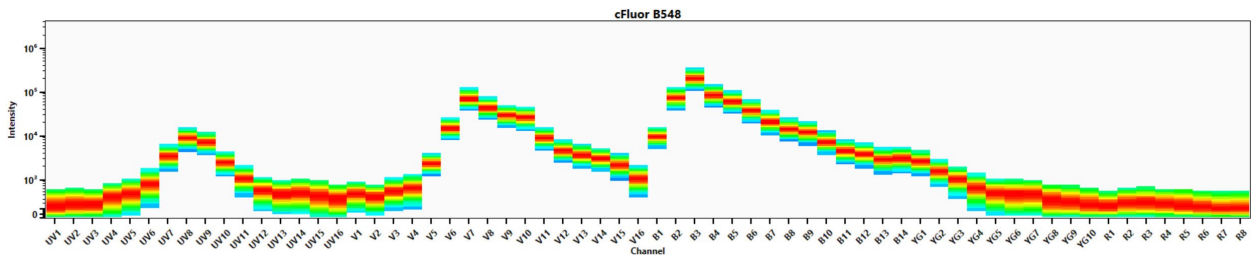

Spectral signature of cFluor® B548 from a Cytek® Aurora 5 laser system equipped with 355, 405, 488, 561 and 640 nm lasers using CytekAssaySetting.

REFERENCES

1. Terry LA, et al. 1990. Immunology. 64:331
2. Hermiston ML, et al. 2003. Annu Rev Immunol. 21:107

**For Research Use Only. Not intended for use in diagnostic procedures.**

\*cFluor® B548 is equivalent to CF® 514 manufactured and provided by Biotium, Inc. under an Agreement between Biotium and Cytek (LICENSEE). The manufacture, use, sale, offer for sale, or import of the product is covered by one or more of the patents or pending applications owned or licensed by Biotium. The purchase of this product includes a limited, non-transferable immunity from suit under the foregoing patent claims for using only this amount of product for the purchaser's own internal research. No right under any other patent claim, no right to perform any patented method, and no right to perform commercial services of any kind, including without limitation reporting the results of purchaser's activities for a fee or other commercial consideration, is conveyed expressly, by implication, or by estoppel.