

cFluor[®] B515 Anti-Human CD7 (CD7-6B7)

PRODUCT DETAILS	
Catalog Number:	R7-20295 (100 tests) R7-20296 (25 tests)
Reactivity:	Human
Clone:	CD7-6B7
Format:	cFluor [®] B515*
Isotype:	Mouse IgG2a, κ
Test Dilution:	5 μL / test
Application:	Flow cytometry
Formulation:	Phosphate-buffered saline, pH 7.2, containing 0.09% sodium azide and 0.2% BSA (BSA Country of Origin USA)
Storage:	2-8°C and protected from light. Do not freeze

PRODUCT DESCRIPTION

The CD7-6B7 monoclonal antibody binds to human CD7, a 40-kDa type I transmembrane glycoprotein in the immunoglobulin superfamily. CD7, also known as GP40, is found on T cells, NK cells, thymocytes, and pre-B cells¹. CD7 is called T cell leukemia antigen for its expression on acute lymphocytic leukemia and some acute myeloid leukemia². Studies show that cross-linking of CD7 significantly increases CD7-associated PI3-kinase activity, suggesting the antigen's functional role in inducing calcium flux in T lymphocytes, and T and NK cell activation³. The antibody was conjugated to a fluorophore and purified by chromatography.


RECOMMENDED USAGE

Each lot of this antibody is quality control tested using flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 μL per 1 million cells in a staining volume of 100 μL. If whole blood is analyzed, then use 5 μL per 100 μL. It is recommended that users titrate the antibody to obtain the optimal result for their specific application.


Please briefly centrifuge the reagent vial before use.

Use appropriate personal protective equipment per the product safety data sheet when using this product.

PRODUCT DATA

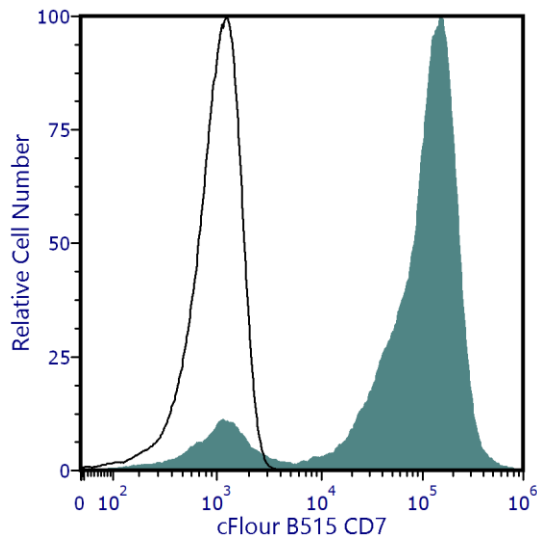

Human peripheral blood stained with cFluor® B515 CD7 (clone CD7-6B7) (filled histogram) or mouse cFluor® B515 IgG2a, κ isotype control (open histogram). Data shown is gated on lymphocytes.

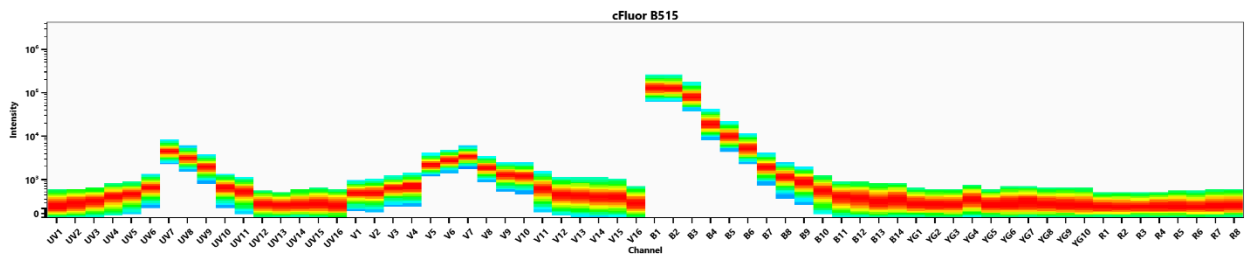

Spectral signature of cFluor® B515 from a Cytek® Aurora 5 laser system equipped with 355, 405, 488, 561 and 640 nm lasers using CytekAssaySetting.

REFERENCES

1. Rabinowich H, 1994. J Immunol. 152(2):517-526.
2. Aruffo A, et al. 1987. EMBO J. 6(11):3313-3316.
3. Lee DM, et al. 1996. Int Immunol. 8(8):1195-1203.

For Research Use Only. Not intended for use in diagnostic procedures.

*cFluor® B515 is equivalent to CF®488A manufactured and provided by Biotium, Inc. under an Agreement between Biotium and Cytek (LICENSEE). The manufacture, use, sale, offer for sale, or import of the product is covered by one or more of the patents or pending applications owned or licensed by Biotium. The purchase of this product includes a limited, non-transferable immunity from suit under the foregoing patent claims for using only this amount of product for the purchaser's own internal research. No right under any other patent claim, no right to perform any patented method, and no right to perform commercial services of any kind, including without limitation reporting the results of purchaser's activities for a fee or other commercial consideration, is conveyed expressly, by implication, or by estoppel.