

cFluor[®] R720 Anti-Human CD16 (3G8)

PRODUCT DETAILS	
Catalog Number:	R7-20005 (100 tests) R7-20006 (25 tests)
Reactivity:	Human, Baboon, Capuchin Monkey, Chimpanzee, Common Marmoset, Cynomolgus, Pigtailed Macaque, Rhesus, Sooty Mangabey, Squirrel Monkey
Clone:	3G8
Format:	cFluor [®] R720*
Isotype:	Mouse IgG1, κ
Test Dilution:	5 μL / test
Application:	Flow cytometry
Formulation:	Phosphate-buffered saline, pH 7.2, containing 0.09% sodium azide and 0.2% BSA (BSA Country of Origin USA)
Storage:	2-8°C and protected from light. Do not freeze

PRODUCT DESCRIPTION

The 3G8 monoclonal antibody binds to IgG receptor III (FcγRIII) that are in two forms: CD16a ((FcγRIIIA) and CD16b ((FcγRIIIB). With 95% sequence similarity, they are a conventional 50-65 kD polypeptide-anchored transmembrane protein and a 48 kD GPI-anchored protein, respectively. CD16a is expressed on NK cells and macrophages while CD16b is expressed on neutrophils^{1,2}. CD16a also plays a crucial role for antibody-dependent cellular cytotoxicity (ADCC) by NK cells³. The antibody was conjugated to a fluorophore and purified by affinity chromatography.


RECOMMENDED USAGE

Each lot of this antibody is quality control tested using flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 μL per 1 million cells in a staining volume of 100 μL. If whole blood is analyzed, then use 5 μL per 100 μL. It is recommended that users titrate the antibody to obtain the optimal result for their specific application.


Please briefly centrifuge the reagent vial before use.

Use appropriate personal protective equipment per the product safety data sheet when using this product.

PRODUCT DATA

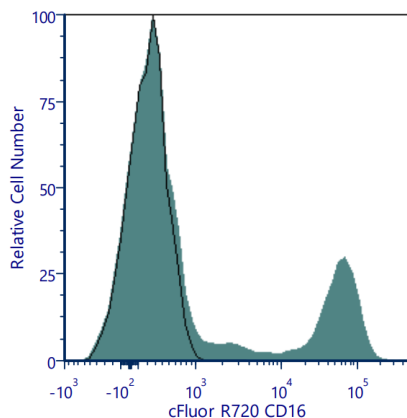

Human peripheral lymphocytes stained with cFluor® R720 CD16 (clone 3G8) (filled histogram) or mouse cFluor® R720 IgG1, κ isotype control (open histogram).

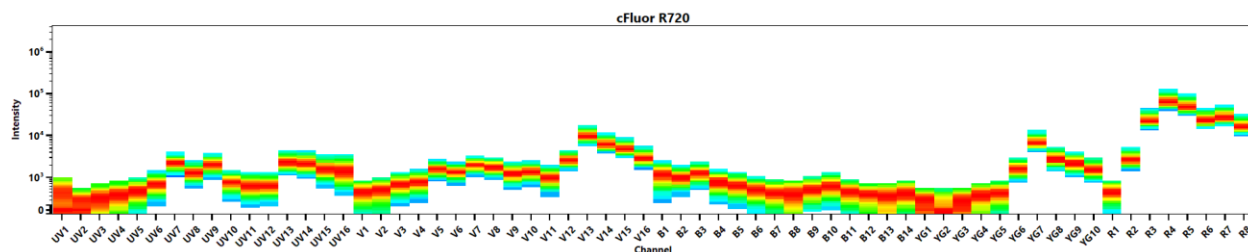

Spectral signature of cFluor® R720 from a Cytex® Aurora 5 laser system equipped with 355, 405, 488, 561 and 640 nm lasers using CytexAssaySetting.

REFERENCES

1. Wirthmueller U, et al. 1992. J Exp Med. 175:1381
2. Smed-Sørensen A, et al. 2008. Blood. 111:5037
3. Wei H Y, et al, 2016. Sci Rep. 6:34310

*cFluor® R720 is equivalent to CF® 700 manufactured and provided by Biotium, Inc. under an Agreement between Biotium and Cytex (LICENSEE). The manufacture, use, sale, offer for sale, or import of the product is covered by one or more of the patents or pending applications owned or licensed by Biotium. The purchase of this product includes a limited, non-transferable immunity from suit under the foregoing patent claims for using only this amount of product for the purchaser’s own internal research. No right under any other patent claim, no right to perform any patented method, and no right to perform commercial services of any kind, including without limitation reporting the results of purchaser’s activities for a fee or other commercial consideration, is conveyed expressly, by implication, or by estoppel.

For Research Use Only. Not intended for use in diagnostic procedures.