

cFluor[®] R720 Anti-Human CD123 (6H6)

PRODUCT DETAILS	
Catalog Number:	R7-20013 (100 tests) R7-20014 (25 tests)
Reactivity:	Human, Rhesus
Clone:	6H6
Format:	cFluor [®] R720*
Isotype:	Mouse IgG1, κ
Test Dilution:	5 μ L / test
Application:	Flow cytometry
Formulation:	Phosphate-buffered saline, pH 7.2, containing 0.09% sodium azide and 0.2% BSA (BSA Country of Origin USA)
Storage:	2-8°C and protected from light. Do not freeze

PRODUCT DESCRIPTION

The 6H6 monoclonal antibody binds to human CD123, a 70-kDa type I transmembrane glycoprotein in the immunoglobulin superfamily, also known as the Interleukin-3 receptor alpha chain. CD123 alone has low affinity to IL-3 and develops high affinity when associated with β c (common chain)¹. The CD123 molecule is expressed predominantly on hematopoietic progenitor cells, macrophages, dendritic cells, basophils, eosinophils, monocytes, megakaryocytes and some B cells^{2,3}. It is also present in myeloid, lymphoid, and vascular endothelial cells⁴. The antibody was conjugated to a fluorophore and purified by affinity chromatography.


RECOMMENDED USAGE

Each lot of this antibody is quality control tested using flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 μ L per 1 million cells in a staining volume of 100 μ L. If whole blood is analyzed, then use 5 μ L per 100 μ L. It is recommended that users titrate the antibody to obtain the optimal result for their specific application.


Please briefly centrifuge the reagent vial before use.

Use appropriate personal protective equipment per the product safety data sheet when using this product.

PRODUCT DATA

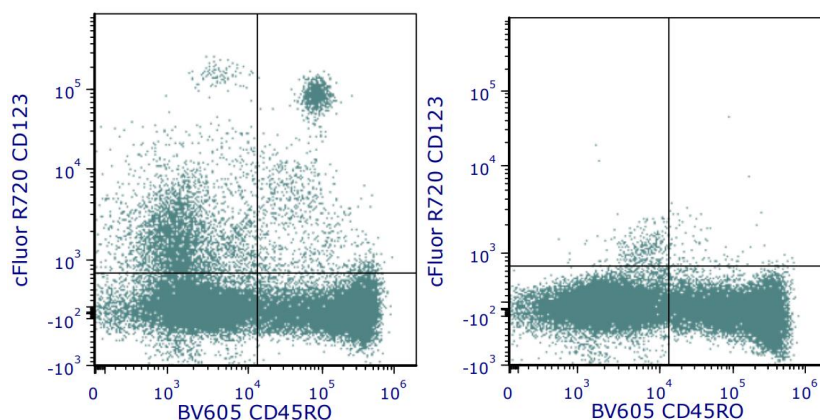

Human peripheral lymphocytes stained with Brilliant Violet™ 605 CD45RO and cFluor® R720 CD123 (clone 6H6) (left) or mouse cFluor® R720 IgG1, κ isotype control (right).

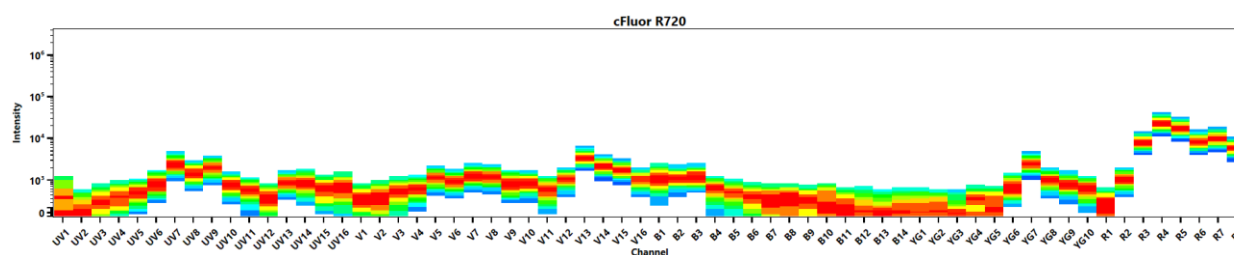

Spectral signature of cFluor® R720 from a Cytek® Aurora 5 laser system equipped with 355, 405, 488, 561 and 640 nm lasers using CytekAssaySetting.

REFERENCES

1. Broughton SE, et al. 2014. Cell Rep. 8:410
2. Macardle PJ, et al. 1996. Cell Immunol. 168:59
3. Yamada T, et al. 1998. J Allergy Clin Immunol. 101:677
4. Sun Q, et al. 1996. Blood. 87:83

*cFluor® R720 is equivalent to CF® 700 manufactured and provided by Biotium, Inc. under an Agreement between Biotium and Cytek (LICENSEE). The manufacture, use, sale, offer for sale, or import of the product is covered by one or more of the patents or pending applications owned or licensed by Biotium. The purchase of this product includes a limited, non-transferable immunity from suit under the foregoing patent claims for using only this amount of product for the purchaser's own internal research. No right under any other patent claim, no right to perform any patented method, and no right to perform commercial services of any kind, including without limitation reporting the results of purchaser's activities for a fee or other commercial consideration, is conveyed expressly, by implication, or by estoppel.

For Research Use Only. Not intended for use in diagnostic procedures.